



Specification

Product Name: Skin Scrubber
Brand: Anymi
Model: PM-2701
Weight: 170g
Power: AC 100-240V, 60Hz

Caution

Never use this device:

- If you are under the age of 12.
- If you have any form of dermatitis.
- If you are prone to allergic reactions or rashes when using skincare products.
- If you have hemophilia.
- If you have liver disorders or any skin disorders related to liver disease or failure.
- If you have pacemaker or any other kind of implanted electronic medical device.
- If you have any metal plates or implants including those in the hips or legs.
- If you are pregnant or breast feeding.
- If you have had or currently have heart problems/heart disease

ProMax Co., Ltd.
ILW00 B/D, 21, Saneop-ro 7 beon-gil, Ojeong-gu, Bucheon-si, Gyeonggi-do, Korea
HQ office: +82-2-839-2011
Website: www.pmkor.com
e-mail: promax@pmkor.com



ANYMI
Respect for human life and Pursuit of beauty

1. COMPONENTS



Device




Charger



Manual

2. NAME OF EACH PART

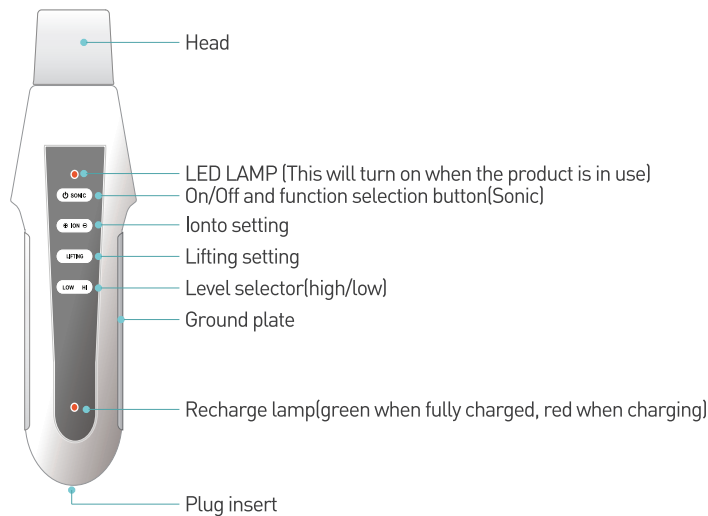

3. FUNCTION

1. Sonic: Removal of dead skin by ultrasonic waves
2. Negative ion: Lets useful ingredient of cosmetics absorbed deep into skin (Iontophoresis)
3. Positive ion: Have skin tone clear by extracting foreign substance in pore (ion discharge/cleansing)
4. Lifting (Works with sonic at the same time): Improves elasticity of skin
5. Low/High: Selects strength of output energy as Low or High
6. Automatic power OFF: Power is automatically switched OFF if it does not operate for one minute after power is switched ON.
7. Safety system: Power is switched OFF automatically when it operates for 12 minutes after power ON. (Prevents overuse)

4. PURPOSE OF USE

1. Removes dead skin by ultrasonic waves (Face, body), however callus like dead skin of heel is excluded.
2. Cleansing by extracting foreign substance in pore (like face)
3. Lets useful ingredient of cosmetics absorbed deep into skin
4. Improves elasticity of skin


5. HOW TO USE

Before use, wash the face cleanly and make it moist. (Do not use it immediately to dry skin)
After cleaning face, please massage skin and open pore, using warm wet towel or steamed towel

1. Ultrasonic waves: Removes dead skin (Use angular side of head corner)
 - 1) Use it after washing the face cleanly using cleansing cosmetics.
 - 2) Press Sonic Button once to turn on device.
 - 3) Select level (strong/weak)
 - 4) Place head of device on skin while thumb and fingers hold the surface of device as most widely as possible.
 - 5) Place 90 degree angular corner of edge of device head on skin and adjust angle of inclination between surface of skin and device as 45 degree. If device is placed on the skin correctly, you can hear and feel the vibration of device while red LED lamp is switched ON.


It means that device operates and it is completely ready for use.
- 6) Move device slowly along with part of body where you want to use device while touching skin to the extent that head of device would not be separated from skin.
- 7) Keep on pushing Sonic Button long to turn off device.
- 8) After skinning, wash the skin and clean head of device with alcohol

Vibration 40,000 times sec

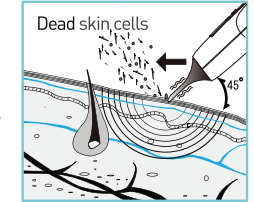



2. Positive ion massage: After washing face, cleansing foreign substance in pore (Cleansing)

- 1) Apply cosmetics for cleansing use or purified clean water properly to the area where you want.
- 2) Press Sonic Button once to turn on device.
- 3) If pressing Ion button once, then positive (+) ion is set.
- 4) Select level (strong/ weak)
- 5) Place head of device on skin while thumb and fingers hold the surface of device as most widely as possible.
- 6) After placing C cut side of edge of device head on skin and adjusting angle of inclination between skin and device as 20 degree, pull it slowly from inside of face (taking nose as a center) to outside. If device is placed correctly, red LED lamp is switched OFF.
- 7) Move it slowly along with area where you want. Device shall be switched OFF automatically if it is left without operation for one minute.
- 8) Press Sonic Button to turn off device.
- 9) Clean part of body with tissue or cotton pad.

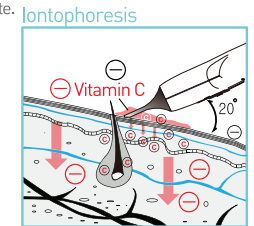

3. Negative ion massage: It helps useful ingredient of cosmetics and Vitamin C to be absorbed deep into skin (Iontophoresis)

- 1) After washing face (inclusive of cleansing), apply cosmetics or Vitamin C properly to the area where you want.
- 2) Press Sonic Button once to turn on device.
- 3) If pressing Ion button once, then negative (-) ion is set.
- 4) Select level. (strong/ weak)
- 5) Place head of device on skin while thumb and fingers hold the surface of device as most widely as possible.
- 6) After placing C cut side of edge of device head on skin and adjusting angle of inclination between skin and device as 20 degree, pull it slowly from inside of face (taking nose as a center) to outside. If device is placed correctly, red LED lamp is switched ON.
- 7) If holding on device in one area for 10 seconds, cosmetics or Vitamin C penetrates into skin. Use device repeatedly on the area where you want. Device shall be switched OFF automatically if it is left without operation for one minute.
- 8) Press Sonic Button to turn off device.

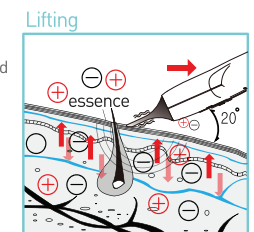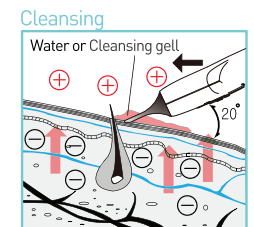

4. Ultrasonic waves and lifting (Repeated operation of negative and positive ion) Concurrent operation: Function that restores elasticity of skin

- 1) Apply nourishing cream (elasticity) or essence properly to the area where you want.
- 2) Press Sonic Button once to turn on device.
- 3) If pressing lifting button, then sonic button is switched ON automatically.
- 4) Select level. (strong/ weak)
- 5) Place head of device on skin while thumb and fingers hold the surface of device as most widely as possible.
- 6) After placing C cut side of edge of device head on skin and adjusting angle of inclination between skin and device as 20 degree, pull it slowly from inside of face (taking nose as a center) to outside. If device is placed correctly, red LED lamp is switched ON.
- 7) Move it slowly along with area where you want. Device shall be switched OFF automatically if it is left without operation for one minute.
- 8) Press Sonic Button to turn off device.